
# SPORTZ-VIBE

MASSAGE THERAPY RUG

*for horses*

**User Manual**

Manuel de l'Utilisateur | Gebrauchsanweisung

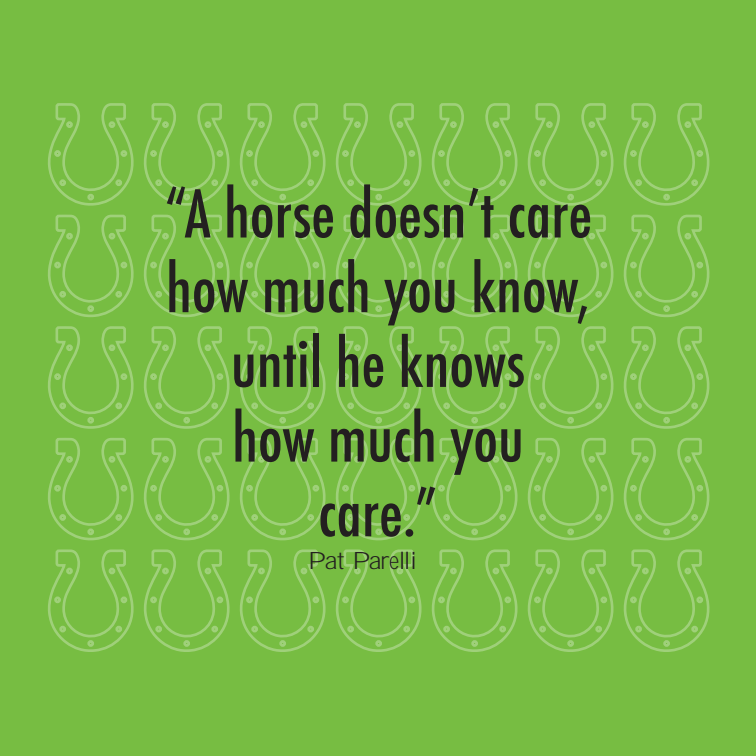**"A horse doesn't care  
how much you know,  
until he knows  
how much you  
care."**

Pat Parelli

## Sportz-vibe for horses | Everything you need to know

If your horse has physical problems that might be aggravated by massage or vibration, consult your veterinarian before using Sportz-Vibe. In some cases, your horse may have an underlying condition which means you shouldn't use Sportz-Vibe, for example in the case of fractures, cancer and certain skin conditions massage may not be an appropriate treatment.

### What comes with the Sportz-Vibe massage rug?

4 panels

1 rechargeable battery

4 panel cables

Battery charger

Lightweight breathable sheet and neck piece

*Extra batteries and panels are sold separately and can be bought from your local stockist.*

## How does it work?

**Warm up:** The increase in blood flow helps to make muscles more elastic and therefore assists you to warm up your horse reducing the risk of injury.

**Recovery:** Massage can help to reduce inflammation by stimulating the lymphatic system to remove toxins. An increase in blood flow also helps bring oxygen and nutrients to the area. This can assist with the healing process.

Your Sportz-Vibe rug has eight different panel positions:

1. Either side of the hind quarters
2. Either side of the saddle area
3. Either side of the shoulders
4. Either side of the neck

There are 3 separate battery positions.

Your rug has four panels and a battery supplied with it. You can move the panels and battery around to treat different areas on your horse.

This allows you to treat hind quarters and the saddle area or the saddle area and shoulders at the same time.

The neck piece has a separate battery pocket and can be treated separately.

## How do I use it?

1. Insert the vibrating panels into the pockets on the underside of the rug. The flat side of the panels should be against your horse's back.
2. Plug the single pin end of the connecting cable into the battery on the right side of the battery pack when you hold it with indicator lights facing you and slide the battery pack into the pocket on the inner of the rug.
3. Push the cables through the elastic loops to keep them covered and hold them in the correct position, the red plug ended cables should go into the panels on the right side of the horse and the black plug ended cables should go into the panels on the left side of the horse.
4. Fasten the surcingle strap closures around the horse's stomach. Use the slide on the strap to adjust the length until it fits snugly.

5. The front of your rug is adjustable and elasticated to give your horse optimum fit. Firstly close the Velcro straps over each other to your required length. When closed connect the male/female front clips to secure the fit at the front. For further adjustments slide the slider on the front straps until the slack is removed from the straps.
6. The rug should be snug and comfortable but not too tight.
7. You can access the battery on your horse's rug by opening the zipped pocket on the outside of the rug. To switch the rug on hold down the on/off button on the outside of the rug for 3 seconds.
8. The rug will turn on at setting 1 and will automatically switch off after 20 mins. During setting 1 the light on your battery will be green.
9. To increase the setting from setting 1 to setting 2 press the on/off button once and it will increase to setting 2. Setting 2 is the strongest setting, this will automatically switch off after 20 mins. The button on your battery during setting 2 will be red.
10. There are three vertical lights on the side of your battery, this shows you how much life is left in your battery. Three lights being the strongest, one light being the weakest.

## What setting should I use?

You should always start your horse on setting 1 to introduce them to the rug. Setting one is gentle and the alternations between the panels is 2 seconds, this allows your horse to adapt to the massage rug.

**Setting 1** is an ideal setting for warm up, recovery after exercise and rehabilitation. The massage action can reduce muscle tension allowing blood flow to return to the muscle, helping to reduce lactic acid build-up and stimulate repair.

**Setting 2** is to be used if your horse has no physical issues, it will help to reduce muscle tension and help create heat in the muscles, making it an ideal setting before you ride or as a treatment for general well being.

## Washing your Sportz-Vibe rug


Remove the four panels, cables and battery from your rug and wash on a 30 degree wash. Allow to drip dry. Do not tumble dry or dry clean.

## Storing your Sportz-vibe rug

Keep your rug in the bag provided along with your charger to keep all the components together.
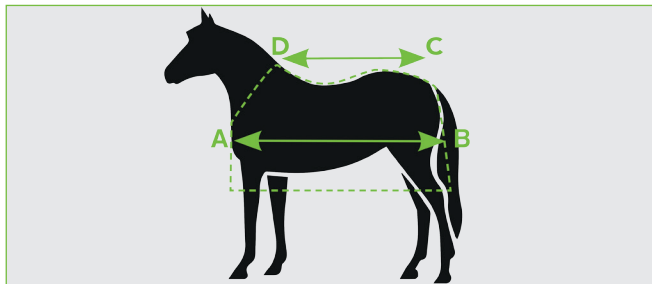
## Rug Measurements

Size range

SIZES		A - B	C - D	INCHES
		ALL DIMENSIONS ARE IN CM		
S		160 - 175	110 - 125	5'3 - 5'9
M		180 - 198	130 - 145	6'0 - 6'6
L		206 - 221	150 - 165	6'9 - 7'3

## User Instructions

Insert battery into battery pocket and seal velcro.

Guide all four cables through elasticated keeper.

Guide one cable of each colour into opposite keepers.

This is how your rug should be laid out.

Once rug is on your horse open battery pocket using outside zip.

Insert panels into pockets.

Panels on right hand side are fitted with red cables.

Panels on left hand side are fitted with black cables.

Hold down power button for 3 seconds until light turns green.

The rug will turn on at setting 1 and will automatically switch off after 20 mins.

## How to charge the battery

Your battery has been partially charged at the factory, but you may need to recharge it before you can switch on your device for the first time.

If not using your battery for more than 4 weeks continue to charge the battery weekly, otherwise it will drain and be unable to recharge.

**USB charging:** your device can be charged from a USB port eg: a laptop, Connect the charger to a wall outlet, connect the charger to the battery, the connection is at the bottom of the battery.

Battery will fully charge in two hours. Do not leave on charge overnight. Disconnect the charger from the wall outlet when the charger is not in use.

Your battery should last for 2 hours once charged.

## Important points to remember

The Sportz-Vibe massage rug is designed to reduce muscle tension. It is not a medical device.

When using your rug make sure your horse has access to plenty of clean drinking water.

To allow your horse to become used to the massage rug, begin by using it once a day on setting 1 for the initial two to three days, then you may increase the sessions per day to a max of 3.

If you wish to change settings, press your battery switch once more.

If you need to replace any parts contact the retailer where you purchased your rug or email us at [sportzvibe@horseware.com](mailto:sportzvibe@horseware.com)

Do not use on fractures as it may create swelling.

Do not use on open wounds as it may spread infection.

Do not use if your horse is unwell and has a temperature.

**Do not use more than 3 times a day.**

If you are unsure please consult your veterinarian.

If have any queries on rug usage please contact us at [sportzvibe@horseware.com](mailto:sportzvibe@horseware.com)

## User guide safety

Only qualified personnel may repair this product.

Keep your electronics dry. They are not water resistant.

Do not store the electronic components in extreme hot or cold temperatures.

Do not drop, knock or shake device.

Do not allow your horse to chew or stand on device as it contains delicate circuit boards.

Use battery and charger for their intended purposes only. Improper use, or use of unapproved batteries or chargers may be harmful and will invalidate warranty.

Never use a damaged charger.

All electrical and electronic products and batteries must be taken to separate collection at the end of their working life. This requirement applies in the European Union.

Do not dispose of these products as unsorted municipal waste.

## Warranty

The warranty period begins from the time of original purchase by the first end user.

The product has a 6 month warranty.

The limited warranty period will not be extended or renewed due to subsequent resale.

This limited warranty does not cover defects caused by rough handling or chewing from the animal, defects caused by sharp items or damage caused by misuse of the product.

This limited warranty does not cover defects or damage caused by the fact that the product was connected to or used with any product not supplied by Horseware®.

This limited warranty is not enforceable if the product has been exposed to dampness, extreme thermal or environmental conditions or to spillage of food, liquid or chemical products.

This limited warranty is not enforceable if the product has been modified or repaired by an unauthorised person.

**If you require any further information or have any questions please contact [sportzvibe@horseware.com](mailto:sportzvibe@horseware.com)**

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